

ImpactIQ Wizard

ImpactIQ Wizard - a software package providing automated evaluation with respect to image quality assurance of CT and CBCT equipment.


ImpactIQ Wizard is a software tool allowing fast and quantitative assessment of image quality in X-ray computed tomography.

The software can be installed on any Windows™ based computer system. Load a specific dicom image set - and only a few clicks are required to receive a printable report with all necessary image quality parameters. The report can be exported easily as pdf.

ImpactIQ Wizard is available for several specific phantoms, e.g.:


- DVT-161 phantom in accordance with DIN 6868-161
- dentalCBCT phantom (BASIC) and (EXPERT)

Softwareservice and -support will be provided by our partner:

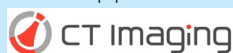


info@ct-imaging.de **www.ct-imaging.de**


*Please note: if using the software, your systems FOV has to completely cover the phantom in height, otherwise a correct automated evaluation is not possible.
Please ask for more information: info@qrm.de*

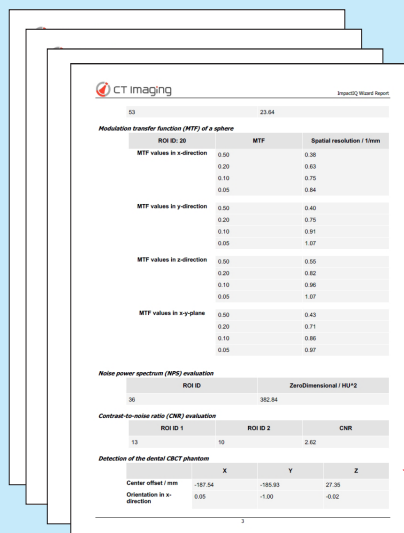


Print report as ,pdf

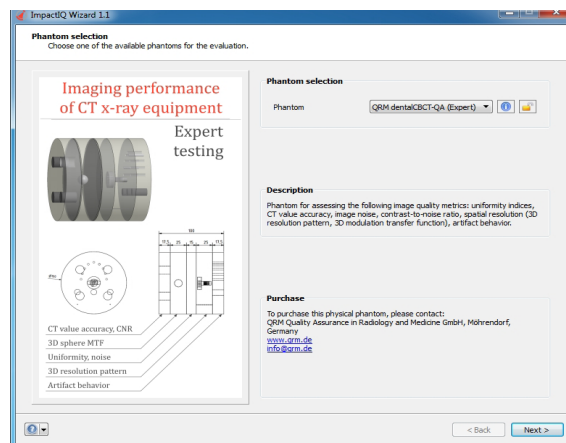

Start ImpactIQ Wizard and select dedicated phantom

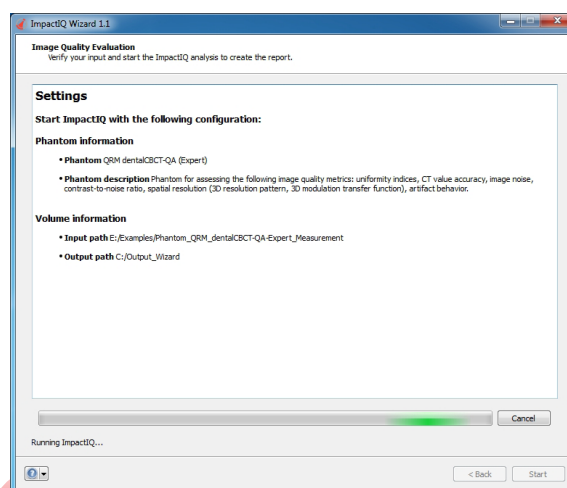

Select DICOM files and start evaluation

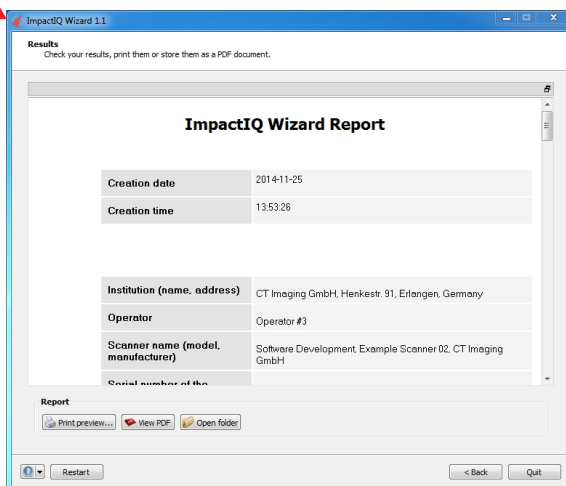

Receive results within seconds