

Pediatric Phantoms

Anthropomorphic thorax and abdomen phantoms in the size of an infant up to a young child.
An excellent tool to validate IQ parameters or dose in a semianthropomorphic surrounding.

QRM-Pediatric-Phantoms providing an anthropomorphic shape of children at different age. These phantoms are giving many opportunities for realistic measurements of dose and image quality in X-ray diagnostics.

The phantoms are designed to provide realistic organs and shapes and offer the physical X-ray attenuation properties of human tissues.

QRM-Pediatric-Thorax comprises artificial lung lopes, a spine insert, and shell of soft tissue equivalent material.

QRM-Pediatric-Abdomen comprises a fixed liver and spleen insert, a spine insert, and shell of soft tissue equivalent material. It is available as well in a blank version (without liver and spleen).

The plastics used in these antropomorphic phantom mimic human tissues with respect to density and X-ray attenuation characteristics of the human thorax and abdomen.


Optional available are a central bore to house an optional insert (mediastinal/stomach region).

For dosimetry purposes 5 bores to cover dose chambers are available as well. The bores are orientated as specified for CTDI measurements.


Ionization chamber adapters are available on request.



1 year, 6 years and young adult Thorax-Phantoms

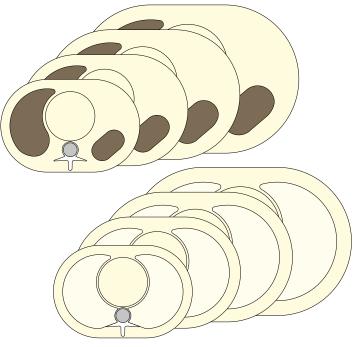


Pediatric phantoms at different age and size

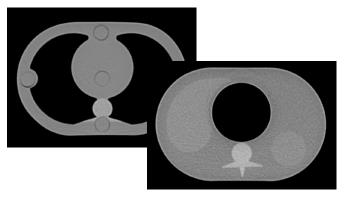


CT-scans of the phantoms



Pediatric Phantoms

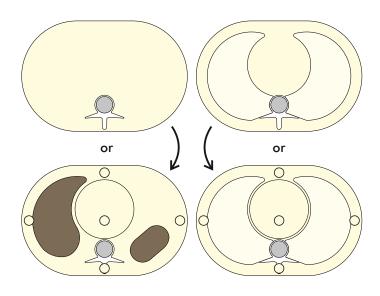

Different versions of phantoms available: standard, additional organs, central exchangeable inserts or dose bores

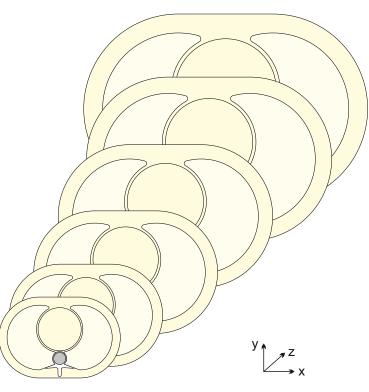

different dimensions available

Specifications

Materials:	tissue	35 HU (at 120 kV)
	lung	app800 HU
	liver / spleen	. 60 HU (at 120 kV)

spine bone like*(

*(offers attenuation properties as trabecular and cortical bone of the specific age

Measures of the phantoms:

Thora Abdor		z [mm]	circumf. [mm]	ВМІ
Newbo	orn 135x90	150	330	
1 year	150x100	150	415	17
3 years	s 180x120	150	497	16
6 years	s 210x140	150	580	16
12 yea	rs 240x160	150	663	18
15 yea	rs 270x180	150	745	20
18 yea adult	300x200	100	828	21

Basic version of Pediatric Thorax:

• lung lobes, spine

Optional:

- dose bores
- central insert (blank)

Basic version of Pediatric Abdomen

• spine

Optional:

- organs (liver/spleen)
- dose bores
- central insert (blank)